
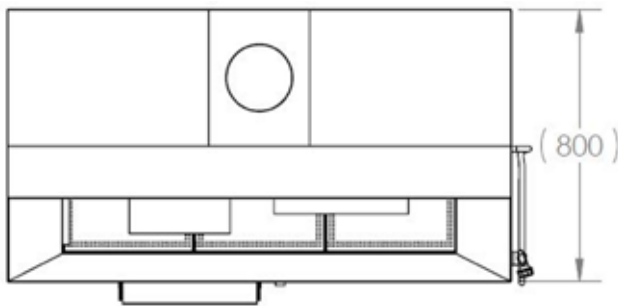




APZEM ELEVATING DUAL DRAFT GROSSING STATION

(DDE Series)

APZEM Elevating series Grossing Station were designed to provide Sturdy, customizable and affordable workstation for medium to large Labs. The Elite Elevating series workstation comes with higher suction capacity, ergonomically designed and additional height adjustable features that assists in high productivity grossing as well as protects you from inhaling harmful Formalin fumes.

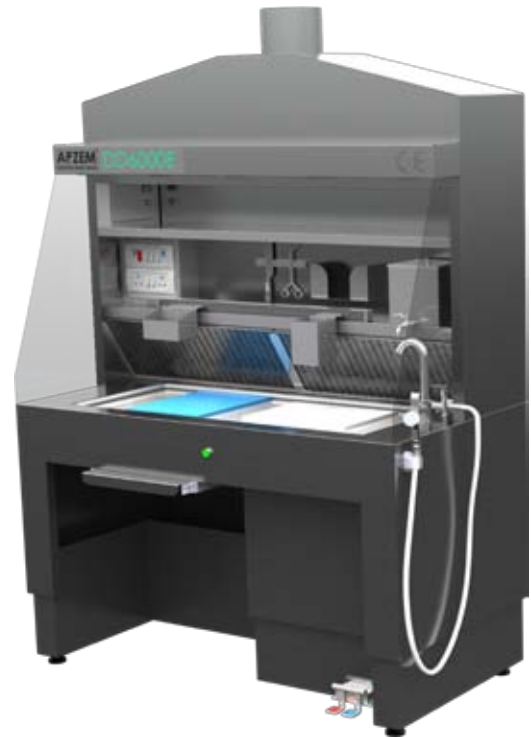

MODEL NO	SIZE (LENGTH)
DD5000E	5 FEET
DD6000E	6 FEET

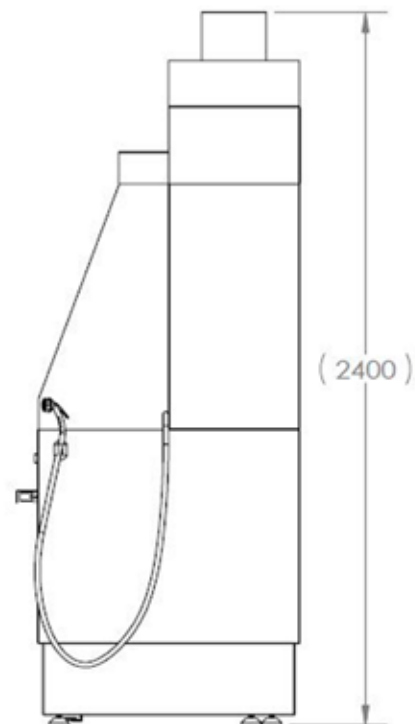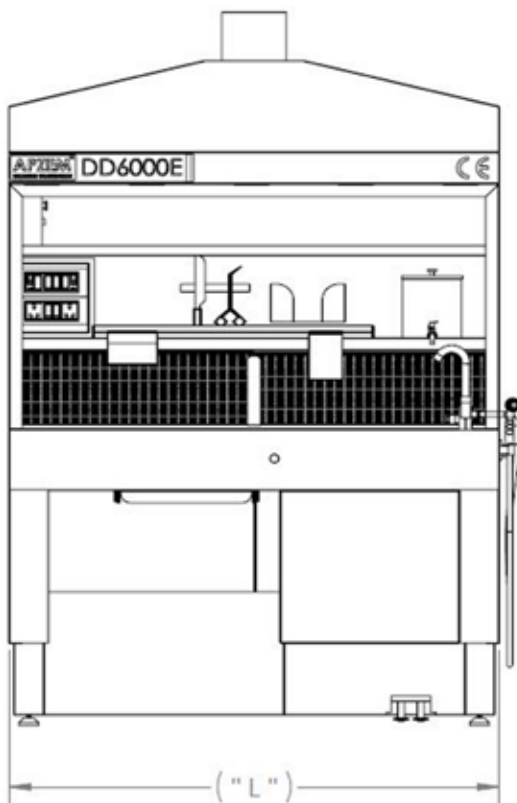

STANDARD DESIGN FEATURES

- ▶ All Stainless steel Construction made of 16 Gauge SS 304.
- ▶ Removable perforated worktop.
- ▶ Inbuilt ventilation system capable of removing all formalin fumes from work top and from Technician's breathing zone.
- ▶ Deep Sink with Cold water mixer Faucet.
- ▶ Back Draft fume Removal system.
- ▶ Powerful LED Lighting.
- ▶ Perimeter rinsing for easy cleaning of work surface.
- ▶ C-Fold paper towel Holder.
- ▶ Polyethylene Dissecting Board for Dissecting and for Photography.
- ▶ Magnetic Tool holder.
- ▶ Cleanable Pre Filter to remove fine particles.
- ▶ Worktable Height : 800mm (normal position).
- ▶ Height can be adjusted from 800mm to 1050mm electro hydraulically with feather touch control switch.
- ▶ Additional shelf for storing.
- ▶ Hand rinsing assembly.
- ▶ Perimeter Rinse assembly.
- ▶ High capacity 0.55 HP waste disposal unit.
- ▶ Transparent Lexan air deflector panels.
- ▶ 40mm height adjustable leg for accurate leveling of work surface.
- ▶ Formalin Dispensing system.

OPTIONAL FEATURES AND ACCESSORIES

- ▶ Magnifying Lamp fixture for mounting 3 axis adjustable magnifying Lamp.
- ▶ Digital camera holder.
- ▶ Safety splash shield.
- ▶ Foot operated mixer faucet.
- ▶ Removable Ventilation Trash cabinet with wheel arrangement.
- ▶ Adjustable camera holder arrangement.
- ▶ Adjustable Shelving.
- ▶ Higher grade, SS 316 Construction for better chemical resistance and greater life.
- ▶ Optional Carbon filter / potassium permanganate filter.
- ▶ Pull out writing pad with storage.

** Chose any of 48 accessories listed at end..