
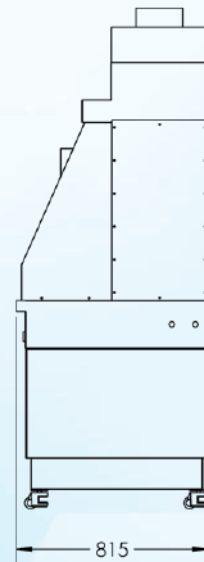
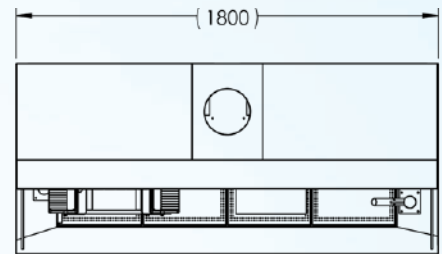


APZEM ELEVATING DUAL DRAFT GROSSING STATION

(DDE SERIES)

APZEM Elevating series Grossing Station were designed to provide Sturdy, customizable & affordable workstation for Elite Labs large Labs. It comes with height adjustable feature, Supports multiuser operation, Inbuilt Gross Imaging System & Formalin Monitor. The Elite Elevating series workstation comes with higher suction capacity, ergonomically designed & additional height adjustable features that assists in high productivity grossing as well as protects you from inhaling harmful Formalin fumes.

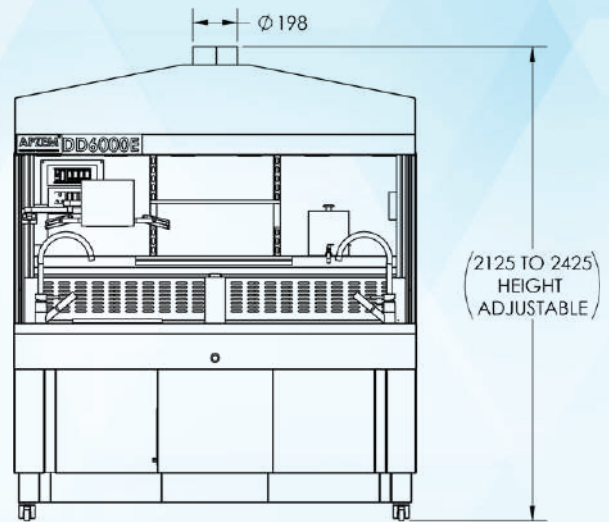

MODEL NO

SIZE (LENGTH)

DD5000E

5 FEET

DD6000E

6 FEET

Specification

Model	DD6000E
MOC	SS304 (Standard) / SS316 (Optional)
Ventilation	Backdraft
Dimension	Minimum : 1850mm Length x 825mm Width x 2150 mm Height Maximum : 1850mm Length x 825mm Width x 2450 mm Height
Work Surface	900mm Minimum height Height can be increased upto 1200mm by feather touch controls
Sink	380mm length x 350mm width x 180mm deep sink
Filters	Water washable polyester prefilters along with optional custom designed filter
Electrical	220 - 240 VAC 50/ 60 Hz, 1 Phase, 110 V supplied on request
Plumbing	1/2" industrial grade 1meter flexible hose supplied for hot and cold water connection 1 1/2" drain hose (waste disposal)
Waste Disposal	0.55 HP High capacity waste disposal unit.
Exhaust air capacity	1800 CMH
Perforated working area	1450mm length x 470mm width
Exhaust Duct Diameter	200mm diameter

Standard Features



Anti Corrosive Construction



Perforated Worktop



Inbuilt Ventilation



Spacious Sink



Backdraft Extraction



Powerful Lights



Boundry Rinse



Towel Holders




Laser Etched Cutting Board




Instrument Holder (Magnetic)




Custom Filters

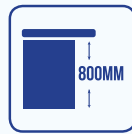

Ergonomic Design

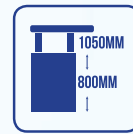

Height Adjustable

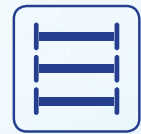

Spacious Shelves



Hand Rinse



Bio Decomposer (Heavy Duty)



Transparent Panels



Work floor Levelling



Formalin Dispenser



Safety Shield



Hands Free Control



Pullout Pad



Formalin Monitoring System



Adjustable Shelves



Path Imaging System

Optional Features



Magnifying Lamp



MOC: SS316



Optional Carbon Filter



Trash Cabinet